


## Product specifications

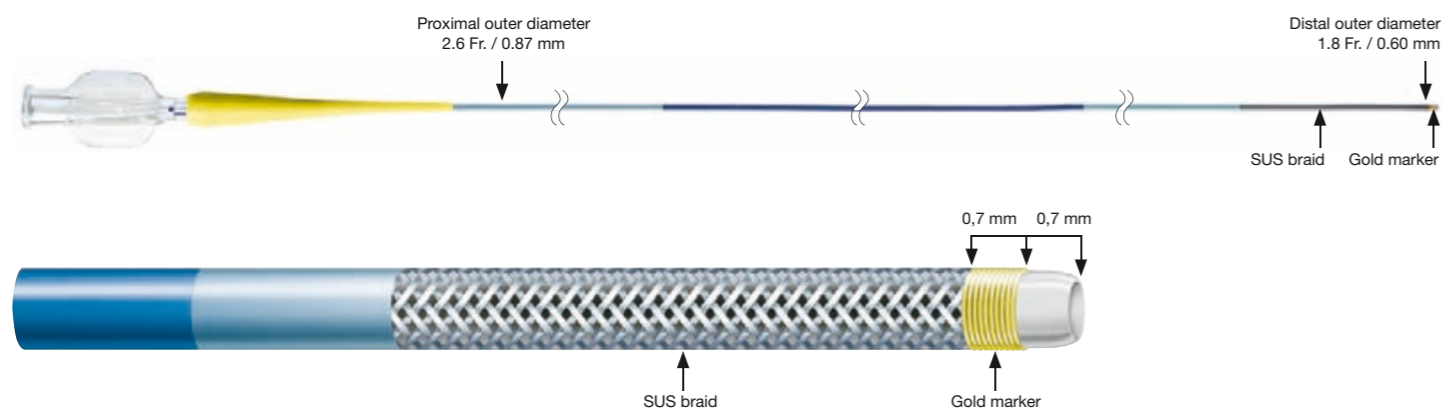

## Coronary Micro-Guide Catheter

## Ordering information

Product code	Distal outer	Proximal outer	Distal inner	Proximal inner	Effective catheter length	Max. guide-wire size
NC-F863A	1.8 Fr. / 0,60 mm	2.6 Fr. / 0,87 mm	0.018" / 0,45 mm	0.021" / 0,55 mm	130 cm	0.014" / 0,36 mm
NC-F865A					150 cm	

Designed for greater  
crossing success  
in complex procedures



**Terumo Corporation**  
44-1, 2-Chome  
Hatagaya, Shibuya-Ku  
Tokyo 151-0072, JAPAN  
Tel: +81/3-3374-8111  
Fax: +81/3-3374-8196  
www.terumo.com

### European Offices

**Terumo Europe N.V.**  
Interleuvenlaan 40  
3001 Leuven  
Belgium  
Tel: +32/16.38.12.11  
Fax: +32/16.40.02.49

**Terumo Europe N.V.**  
Liaison office  
Interventional Systems  
Immeuble Krebs  
8, Rue Jean-Pierre Timbaud  
Cs 10720  
78067 Saint Quentin En Yvelines  
Cedex  
France  
Tel: +33/1 39.30.95.75  
Fax: +33/1 39.30.95.70

**Terumo Europe N.V.**  
Moscow Representative office  
4, Ivana Franko St  
121108 Moscow  
Russia  
Tel: +7/495.234.3303  
Fax: +7/495.234.3303

**Terumo Europe N.V.**  
Benelux Sales Division  
Interleuvenlaan 40  
3001 Leuven  
Belgium  
Tel: +32/16.39.25.80  
Fax: +32/16.39.25.99  
The Netherlands  
Tel: 0800/022.03.96  
Fax: 0800/022.04.14

**Laboratoires Terumo France S.A.**  
Parc Ariane - Bât. Uranus  
1 Rue Hélène Boucher  
78284 Guyancourt Cedex  
France  
Tel: +33/1 30.96.13.00  
Fax: +33/1 30.43.60.85

**Terumo Deutschland GMBH**  
Hauptstrasse 87  
65760 Eschborn  
Germany  
Tel: +49/61.96.80.230  
Fax: +49/61.96.80.23.200

**Terumo Deutschland GMBH**  
Zweigniederlassung Wettingen  
Bahnhofstrasse 88  
Ch-5430 Wettingen  
Switzerland  
Tel: +41/56.419.10.10  
Fax: +41/56.419.10.11

**Terumo Europe España SL**  
Avenida Juan Carlos I, N° 13 - 7ª  
Planta  
Edificio Torre La Garena  
28806 Alcalá De Henares (Madrid)  
Spain  
Tel. General: +34/902.10.12.98  
Cust. Serv.: +34/902.20.31.40  
Fax General: +34/902.10.13.58  
Cust. Serv.: +34/902.10.17.57

**Terumo Italia SRL**  
Via Simone Martini 143/145  
I-00142 Roma  
Italy  
Tel: +39/06.519.61.420  
Fax: +39/06.503.04.07

**Terumo Sweden AB**  
Nya Varvet, Byggnad 90  
Sven Källfelts Gata 18  
Se-426 71 Västra Frölunda  
Sweden  
Tel: +46/31.748.58.80  
Fax: +46/31.748.58.90

**Terumo Denmark**  
Filial Of Terumo Sweden AB  
Tel: +45/7020.93.80  
Fax: +45/7020.94.80

**Terumo UK LTD**  
Tamesis, The Causeway  
Egham  
Surrey  
TW20 9AW  
United Kingdom  
Tel: +44/1784.476.200  
Fax: +44/1784.476.201



ISI39GB1208MDM0510JE

Finecross® MG coronary is the new and innovative micro catheter specifically developed for coronary interventions.

Finecross® MG is designed to provide outstanding guide wire handling and support without compromising the crossability.

The new Finecross® Micro catheter is your best PCI partner to make your work easier in very challenging cases involving, tortuous vessels, CTO, bifurcations, highly stenosed and calcified lesions.


# Guide wire handling

## Innovative technologies for superior\* guide wire handling and optimal wire support

### Superior\* guide wire handling

#### Large inner diameter

Finecross® MG provides an extra large inner lumen  
 Distal : 0.018" (0.45mm)  
 Proximal : 0.021" (0.55mm)  
 This reduces friction and ensures an optimal guide wire handling.

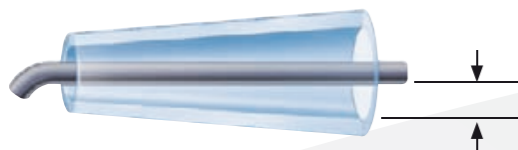


#### Tapered shaft design

Outer distal diameter 1.8Fr (0.60 mm) to proximal outer diameter 2.6Fr (0.87 mm).  
 This provides optimal\* balance between support and crossability.

### Optimal\* wire support

#### SUS braid technology

Finecross® MG employs an innovative stainless steel shaft construction (SUS braid), which covers the entire length of the catheter shaft all the way to the distal tip. This provides strength responsiveness, support, reinforces lumen integrity and improves pushability. This stable platform may help to treat complex lesions including CTO/ Tortuous anatomy.

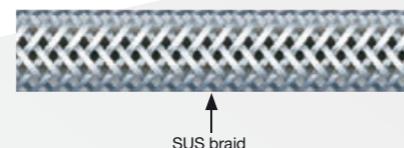

#### PTFE inner layer

PTFE ultra smooth inner surface facilitates manipulation of the guide wires. This provides superior wire handling and facilitates exchange wire during complex procedure.

\* Data on file

# Improved crossability

## Advanced technologies

### 13 cm distal flexible segment



Provides an optimal balance between handling and safety.

### M coat™ hydrophilic coating

The ultra smooth outer surface facilitates advancement through tortuous anatomy and tight lesions.

### Tapered inner and outer diameters

Ensure optimal Guide Wire handling and exceptional crossability.

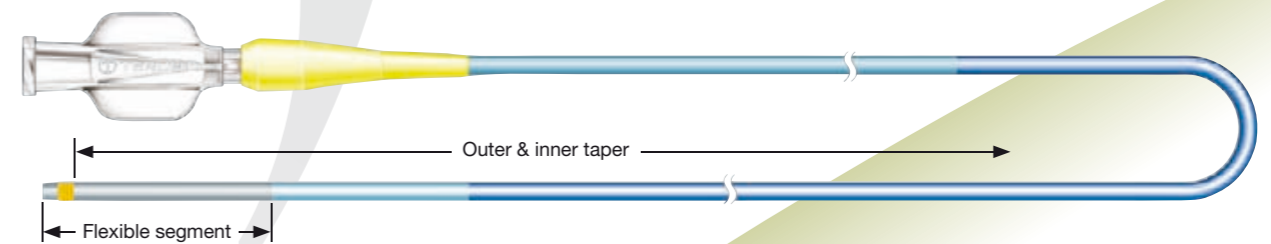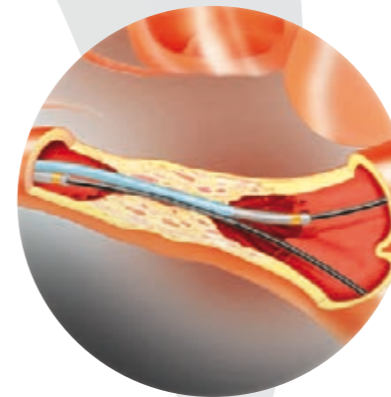

## Reliability you can count on

### Distal gold marker

130 /150 cm lengths

