

# ZIPwire™

Nitinol Hydrophilic Guidewire

Boston  
Scientific

## The Fast Track for Ureteral Access.

*Featuring*  
**Tip and Body Enhancements  
to Stiff ZIPwire Guidewire**



# ZIPwire™ Nitinol Hydrophilic Guidewire


The ZIPwire Guidewire is a kink-resistant nitinol wire with a lubricious hydrophilic coating designed for reliable ureteral access and reduced trauma. The fast track for ureteral access. The ZIPwire Guidewire is available in standard and stiff shaft.

## Lubricious, floppy tip

facilitates reliable access to the ureter and reduced trauma during advancement through the urinary tract

**Durable, covalently bonded hydrophilic coating** designed to navigate ureteral strictures and facilitate tracking of urological instruments

**Excellent torque control**  
when crossing tortuous anatomy

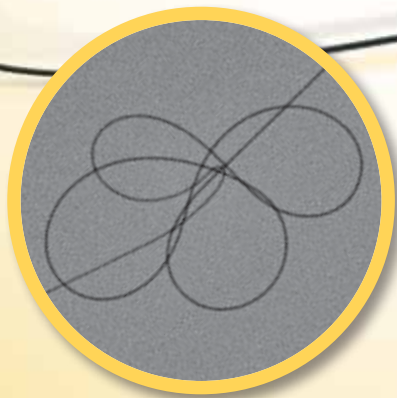

## Highly radiopaque jacket

designed to provide enhanced visibility under fluoroscopy and enable accurate positioning



**Torque vise** may be used for added guidewire steerability and manipulation



**Kink-resistant superelastic nitinol core**  
Designed for ease of advancement


## The Stiffness Option

The Stiff ZIPwire™ Guidewire is designed to increase pushability, help straighten tortuous anatomy and facilitate advancement of devices.

Enhancements have been made to the tip and body of the current Stiff ZIPwire Guidewire to make it stiffer.

	Tip Enhancement	Body Enhancement
<b>.025" Stiff ZIPwire Guidewire</b>	✓	
<b>.035" Stiff ZIPwire Guidewire</b>	✓	✓
<b>.038" Stiff ZIPwire Guidewire</b>	✓	

According to bench testing this table depicts the redesigned segments to be as rigid as compared to other market devices. Hence, our NEW Stiff ZIPwire Guidewires are designed to meet your stiffness needs.

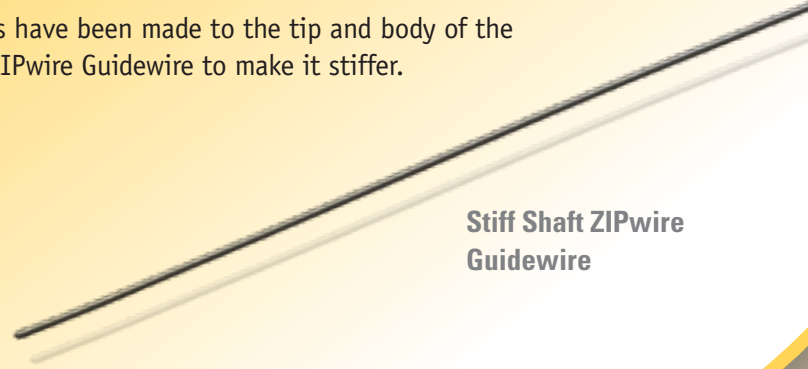


Stiff Shaft ZIPwire Guidewire


## Multiple Tip Styles

Several tip styles are available, including the Bentson-type tip designed to provide the additional tip flexibility for negotiation of narrowed ureters and access beyond severe obstructions.

ZIPwire Guidewire Performance	
<b>Reduced Trauma</b> Soft, flexible tip designed for reduced trauma during advancement through the urinary tract	✓
<b>Lubricious</b> Hydrophilic coating designed for smooth advancement and easy tracking of devices	✓
<b>Kink-Resistant</b> Superelastic nitinol core designed to resist kinking	✓
<b>Steerable</b> Nitinol core and polyurethane jacket provide torqueability for excellent control and positioning	✓
<b>Radiopaque</b> Highly radiopaque jacket designed for enhanced fluoroscopic visibility	✓

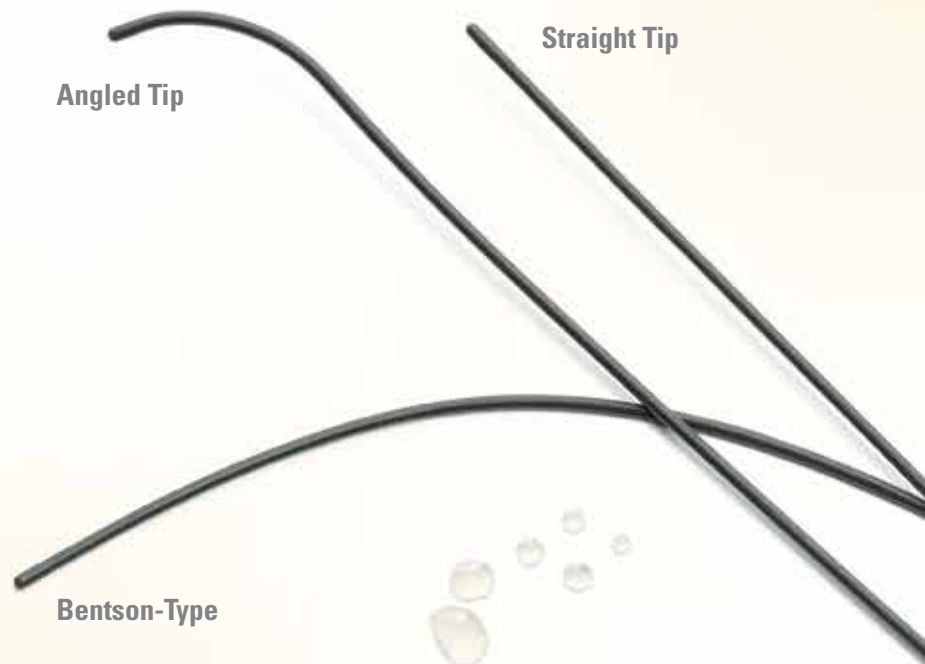

ORDER NUMBER	DIAMETER	LENGTH	TIP	SHAFT STIFFNESS	TAPER	QUANTITY
<b>Standard ZIPwire Guidewire</b>						
M006630200B1	0.018 in	150 cm	Straight	Standard	3 cm	Box 5
M006630201B1	0.025 in	150 cm	Straight	Standard	3 cm	Box 5
M006630202B1	0.025 in	150 cm	Angled	Standard	3 cm	Box 5
M006630212B1	0.032 in	150 cm	Straight	Standard	3 cm	Box 5
M006630213B1	0.032 in	150 cm	Angled	Standard	3 cm	Box 5
M006630205B1	0.035 in	150 cm	Straight	Standard	3 cm	Box 5
M006630206B1	0.035 in	150 cm	Angled	Standard	3 cm	Box 5
M006630208B1	0.038 in	150 cm	Straight	Standard	3 cm	Box 5
M006630209B1	0.038 in	150 cm	Angled	Standard	3 cm	Box 5
M006630210B1	0.038 in	260 cm	Straight	Standard	3 cm	Box 5
<b>Stiff ZIPwire Guidewire</b>						
M006630216B1	0.025 in	150 cm	Straight	Stiff	3 cm	Box 5
M006630217B1	0.025 in	150 cm	Angled	Stiff	3 cm	Box 5
M006630222B1	0.035 in	150 cm	Straight	Stiff	3 cm	Box 5
M006630223B1	0.035 in	150 cm	Angled	Stiff	3 cm	Box 5
M006630225B1	0.038 in	150 cm	Straight	Stiff	3 cm	Box 5
M006630226B1	0.038 in	150 cm	Angled	Stiff	3 cm	Box 5
<b>Bentson-Type ZIPwire Guidewire</b>						
M006630214B1	0.035 in	150 cm	Straight	Standard	Bentson 8 cm	Box 5
M006630203B1	0.035 in	150 cm	Angled	Standard	Bentson 8 cm	Box 5
M006630215B1	0.038 in	150 cm	Straight	Standard	Bentson 8 cm	Box 5
M006630207B1	0.038 in	150 cm	Angled	Standard	Bentson 8 cm	Box 5
<b>Bentson-Type Stiff ZIPwire Guidewire</b>						
M006630224B1	0.035 in	150 cm	Straight	Stiff	Bentson 5 cm	Box 5
M006630221B1	0.035 in	150 cm	Angled	Stiff	Bentson 5 cm	Box 5
M006630227B1	0.038 in	150 cm	Straight	Stiff	Bentson 5 cm	Box 5
M006630228B1	0.038 in	150 cm	Angled	Stiff	Bentson 5 cm	Box 5
<b>Accessories Guidewire</b>						
M0067301300	Torque Vise Steerable Handle					Box 12

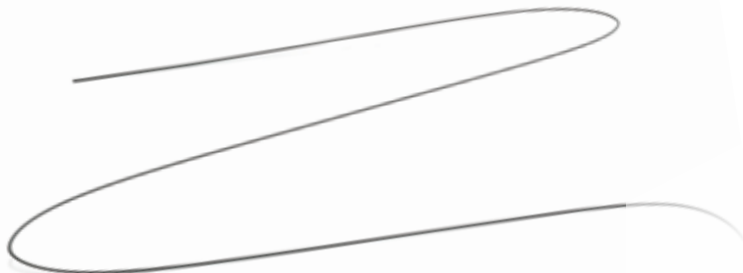

**Boston Scientific**  
Advancing science for life™

Boston Scientific Corporation  
300 Boston Scientific Way  
Marlborough, MA 01752-1234  
www.bostonscientificstone.com

**Ordering Information**  
**1.888.272.1001**

©2014 Boston Scientific Corporation  
or its affiliates. All rights reserved.

URO-135002-AB AUG/2014

All trademarks are the property of their respective owners.

**Caution:** Federal (U.S.) law restricts this device to sale by or on the order of a physician. Please refer to package insert provided with the product for complete Indications for Use, Contraindications, Warnings, Precautions, Adverse Events, and Instructions prior to using this product.

Information for the use only in countries with applicable health authority product registrations.