


Cordis **BRITE TIP**[®] Sheath Introducers

Leading the way
to precision placement

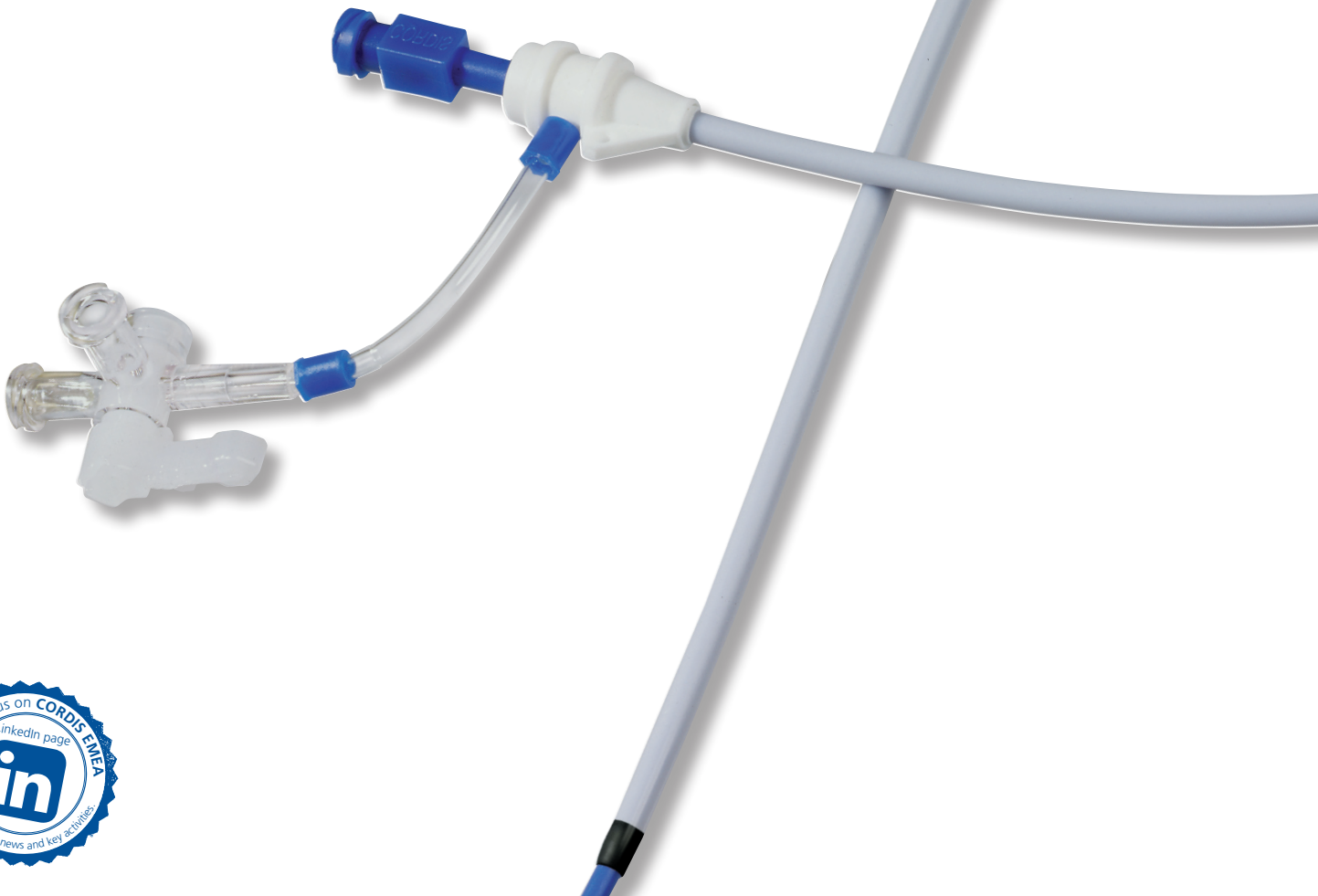

BRITE TIP®

Flexibility and control

BRITE TIP®

Precise, atraumatic placement

Kink resistant cannula

Unique dual layer cannula

Flexibility and support to access tortuous vessels

Thicker cannula wall

Better kink resistance

Silicone coated cannula and vessel dilator


Lubricious, anti-thrombogenic vessel entry



Patented SLIX valve

Six-cut spiral design for hemostatic valve


Free catheter movement without blood reflux or air aspiration

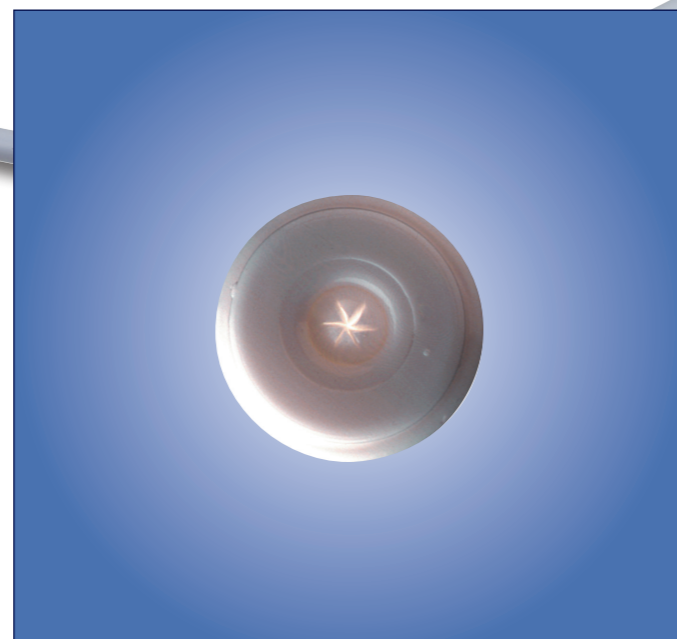

Manicured tip and dilator design

Tapered transitions on sheath and dilator for:

- Effortless vessel entry
- Improved patient comfort
- Potentially reduce groin complications, even when penetrating scar tissue

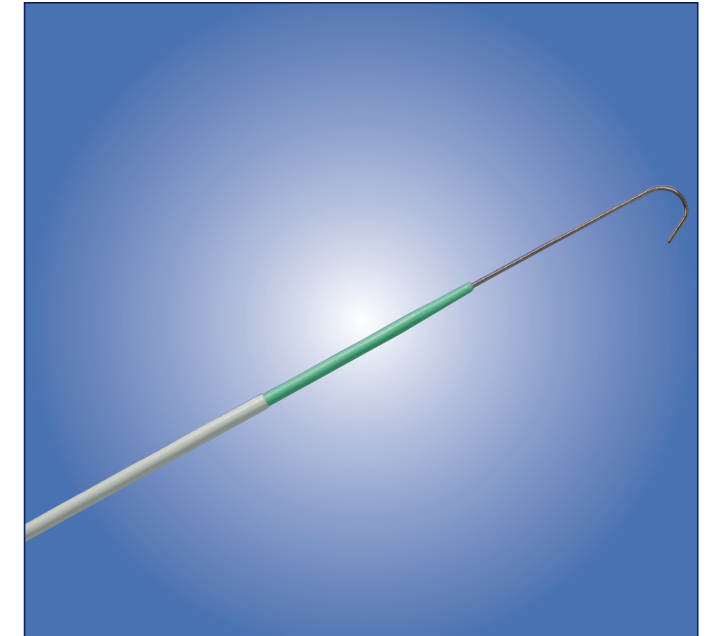

Radiopaque BRITE TIP®

Black tungsten-filled distal tip

Excellent visualization facilitates exact location of sheath's distal tip for accurate interventional positioning



Cordis **BRITE TIP**[®] Sheath Introducers

Leading the way to precision placement

Ordering information

Longer lengths available for a wide variety of interventional procedures



Catalogue Numbers

Cannula Length			5.5cm	11cm	23cm	35cm	45cm	55cm	90cm
French size	Colour code	Dilator Tip (mm)	With mini-wire	With mini-wire	Without mini-wire	Without mini-wire	Without mini-wire	Without mini-wire	Without mini-wire
4	■	25		401-411M	401-423M				
5	■	25		401-511M	401-523M	401-535M	401-545M	401-555M	
6	■	35	401-605M	401-611M	401-623M	401-635M	401-645M	401-655M	401-690M
7	■	35	401-705M	401-711M	401-723M	401-735M	401-745M	401-755M	401-790M
8	■	35	401-805M	401-811M	401-823M	401-835M	401-845M	401-855M	
9	■	35		401-911M	401-923M	401-935M			
10	■	45		401-011M	401-023M				
11	■	45		401-111M	401-123M				

Important Information: Prior to use, refer to the "Instructions for Use" supplied with these devices for indications, contraindications, side effects, suggested procedure, warnings, and precautions. As part of the Cordis policy of continuous product development, we reserve the right to change product specifications without prior notification. For healthcare professionals only.

© Cardinal Health 2016. All right reserved. CORDIS, the Cordis LOGO, BRITE TIP[®] are trademarks or registered trademarks of Cardinal Health.

• EU500-1 10/16

Cordis[®]
A Cardinal Health company